


Drip-Alarm

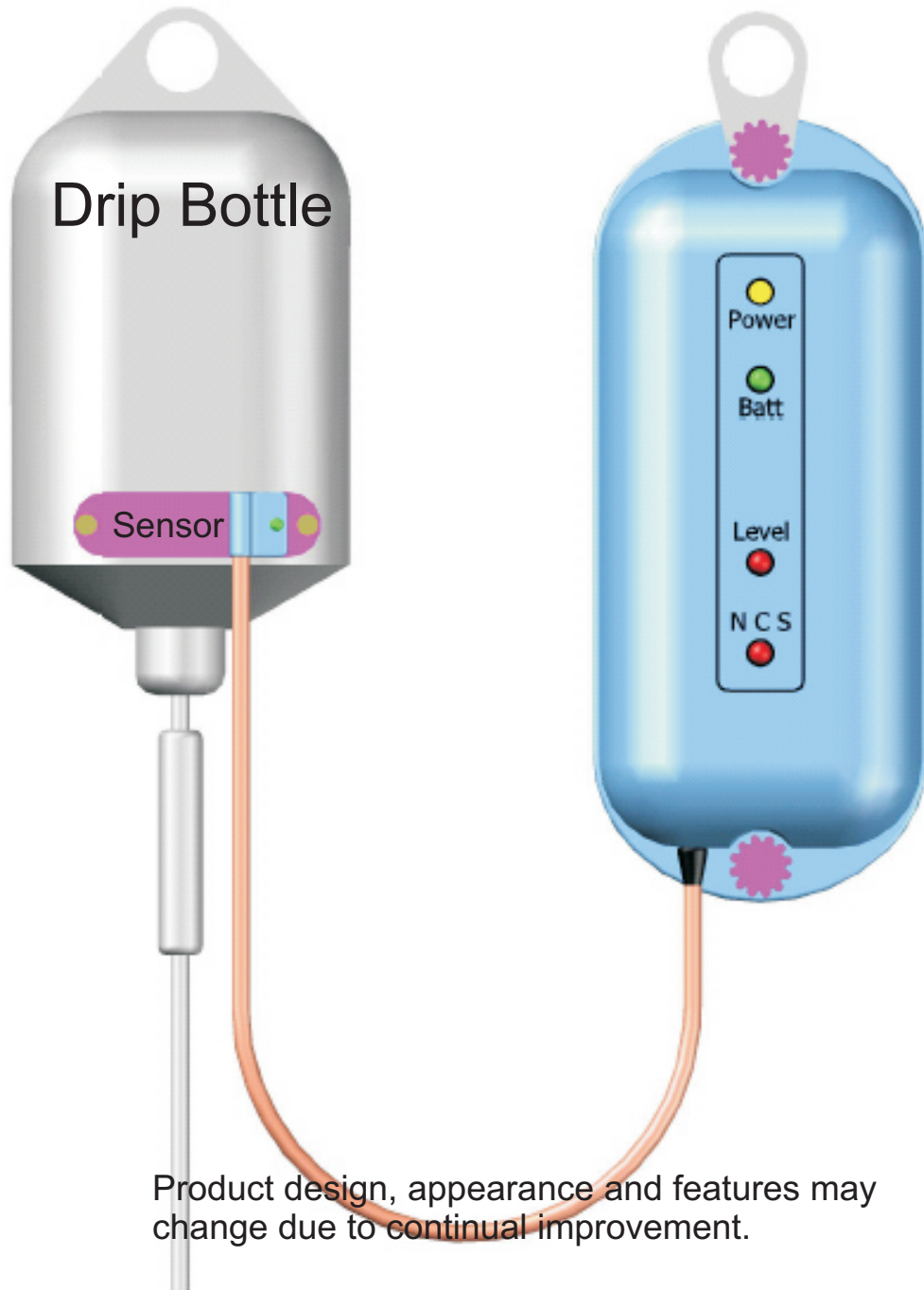

Product design, appearance and features may change due to continual improvement.

The days and nights of keeping vigil on the drip bottle are over. Just fix the Drip-Alarm sensor at the lower part of the drip bottle and relax. Drip-Alarm will start beeping as soon as the liquid falls below the sensing level.

Besides sounding the alert, Drip-Alarm turns on a built-in switch that can be used to operate external alarms or trigger a patient's call on Nightingale type of nurse call systems to remotely alert nursing staff.

Features:

Non-Invasive.

Sensor can be re-used thousands of times.

Sensor can be held in place by a simple rubber band.

Works on plastic or glass bottles.

Easy to use. No controls except On-off switch.

Reassurance light on sensor.

Fool-proof. Alarm sounds if sensor falls off the bottle.

Port for communication to external gadget.

Uses easily available D type torchlight batteries.

Low power consumption ensures long battery life.

Battery level indication by colour changing LED.

Low battery warning by buzzer.

Rugged and elegant ABS plastics body.

Molded ring for hanging next to drip bottle.

Compatible with Nightingale nurse call systems.

Telereach

ISO 9001:2008



#68, 11B Main Road, Jayanagar 5th Block, Bengaluru 560041, India

☎ 80- 26540222, 26649681, 9448494496

telereach.sales@gmail.com